


4D LED LIGHT THERAPY MASK

USER MANUAL

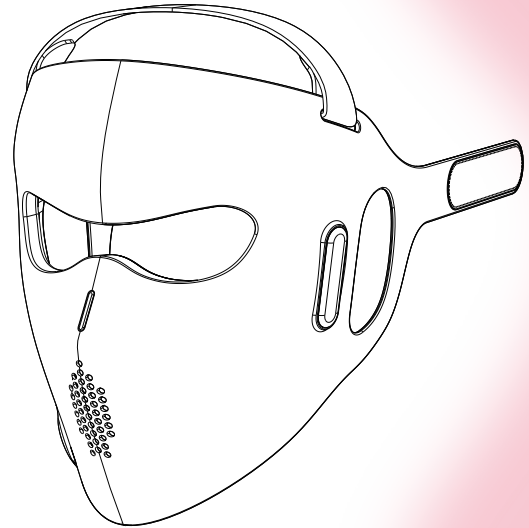


Please read this user manual carefully before beginning and keep it properly for future reference.
Before using this product for the first time, please perform a light sensitivity test.
Please refer to the instructions for details.

Catalogue

Product Description	1
Operating Mode	2
Function	3
LED Arrangement Diagram	4
Product Controller	5
How To Use	6
Product Specifications	6
Installation	7
Safety Cautions	8
Warning	9

Product Description

The 4D LED Light Therapy Mask is an innovative, rechargeable phototherapy beauty device designed for convenient home use. Crafted from medical-grade silicone, this foldable and portable mask features an adjustable buckle system to ensure a comfortable, customized fit for various face shapes. Equipped with 360 high-quality LED chips, the device emits four clinically-proven wavelengths - blue light (415nm), yellow light (590nm), red light (633nm), and near-infrared light (850nm+1072nm) - with four adjustable energy levels through its intelligent controller. Utilizing advanced phototherapy technology, these precisely calibrated light waves penetrate deep into the skin to effectively stimulate collagen and elastin production at the cellular level, delivering visible anti-aging and skin-rejuvenating results. The user-friendly interface offers six preset treatment modes, each with a 10-minute automatic timer for optimal, hassle-free skincare sessions.

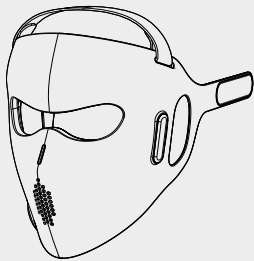

Strap




Controller



Type-C Cable




Carrying Bag

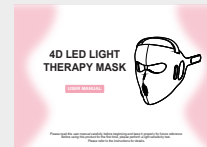

User Manual

Modes

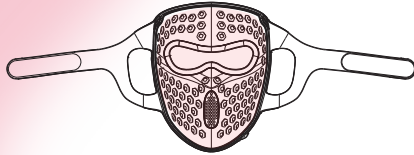

M1
Repairing Mode



M2
Rejuvenation Mode



M3
Anti-Aging Mode



M4
Morning Skincare Mode



M5
Anti-Acne Mode



M6
Bedtime Skincare Mode

Function

6 Modes

M1: Repairing Mode 633nm ■

M2: Rejuvenation Mode 590+633nm ■■

M3: Anti-aging Mode 633+850nm/1072nm ■■

M4: Morning Skincare Mode 590+633+850nm/1072nm ■■■

M5: Anti-Acne Mode 415nm ■

M6: Bedtime Skincare Mode 590+633+850/1072+415nm ■■■■



Yellow Light

Balances skin texture, stimulates red blood cells, reduces redness, swelling, and fine lines



Red Light

Anti-aging, promoting collagen, tightening and increasing elasticity



Near-infrared Light

Repair skin, improve blood circulation, and increase skin moisture




Blue Light

Sterilizes acne producing bacteria, reduces inflammation and scars, and increases skin oxygen

Tips


Before opening the packaging for use, be sure to fully charge the battery of the mask. The charging time is about 3 hours. The upper cover is equipped with an adjustable buckle, allowing you to wear the mask comfortably. To install the mask, thread the strap through the entrance hole on the other side of the mask, wrap the strap around the head and firmly secure it.

LED Arrangement Diagram

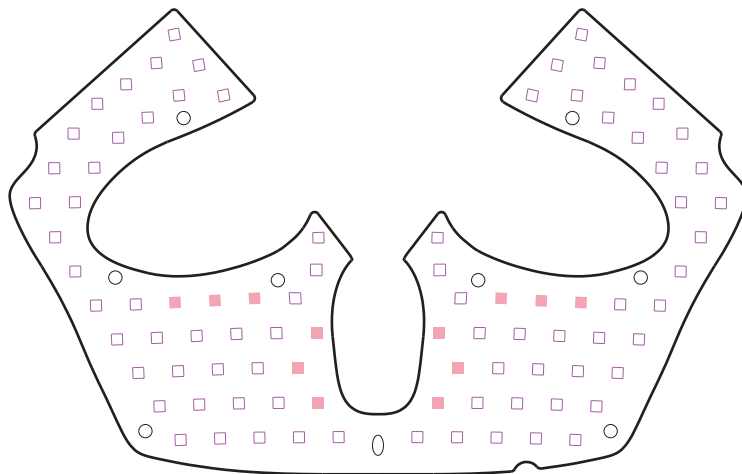

12PCS Pink remark are $415\text{nm} : 590\text{nm} : 633\text{nm} : 850\text{nm} / 1072\text{nm} = 1 : 1 : 1 : 1 / 12\text{pcs}$

Product Controller

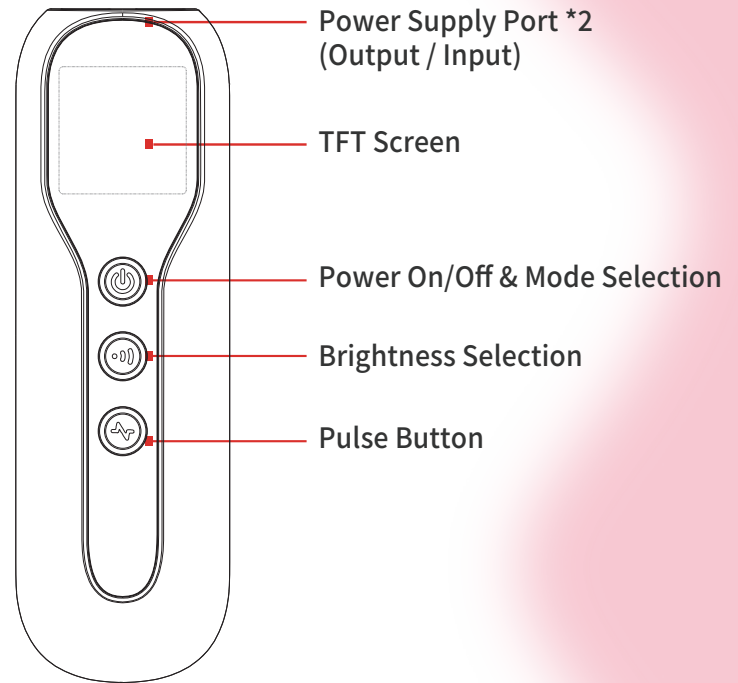

How To Use

- | Power button: Press and hold for about 3 seconds to enter the power on state, and the light will turn on.
- | The brightness selection key: short press to adjust and switch to the 3rd gear energy output mode 50% 75% 100%.
- | Timing: The default timing for each mode is 10 minutes when the mode is selected.
- | Mode selection: Press the power button after the light is turned on. There are 6 modes you can choose.
- | Pulse selection: you can press the pulse function key in any modes.

Product Specifications

Wavelengths	415:590:633:850/1072 nm	LED Color	Blue / Yellow / Red / NIR
LED QTY	90pcs	Power Consumption	G11 6W
Dimensions	345*225*181mm	Charge Part	USB Type-C
Input Voltage	5V 2A/1A	Battery Capacity	4000 mAh (Face piece)
Pulsing	10 Hz	Product Function	6 skin care modes, brightness adjustment, pulse function

Installation


1. Remove dustproof film



2. Put on the adjustable buckle



3. Connect power and controller

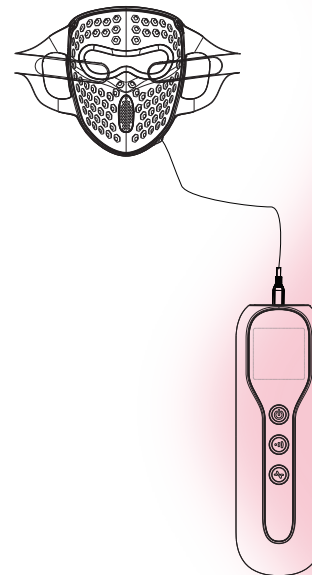

Safety Cautions

Light Sensitivity Test

NOTE: DO NOT USE THIS PRODUCT WITHOUT CONSULTING YOUR HEALTH CARE PROFESSIONAL FOR THE FOLLOWING CONDITIONS

1. People who are pregnant, preparing for pregnancy, or breastfeeding;
2. people with light sensitivity or taking any medication that may cause photosensitivity; people with hereditary diseases of the eye (such as lupus erythematosus, photosensitive eczema or albinism, etc.)

Step 1: Place the device on your arm and turn on the switch to irradiate for 10 minutes. If symptoms such as redness or itching appear on the skin, it indicates that you may have photosensitive skin and it is recommended that you stop using it, or consult a dermatologist before using it again.

Step 2: Before use, please clean your face and keep it dry (without makeup) to prevent photosensitive ingredients that may be contained in cosmetics from affecting the efficacy of skin care.

Step 3: After use, the skin may appear red, which is a normal reaction of the skin, and will return to normal within 1-2 hours. If you feel burning or any other discomfort during the process, please stop using it immediately.

Warnings

1. Please strictly follow the correct method of using this product and do not use the equipment outside of the instructions, otherwise, it may cause machine malfunction or unexpected accidents.
2. Do not use this product to improve any other condition except for indications.
3. Pregnant women, lactating women, and children are not allowed to use this product.
4. If the user suffers from physical, sensory, and mental diseases, or is undergoing dermatology disease improvement, including the onset period (such as lupus erythematosus, photosensitive eczema, or albinism), it is not recommended to use this product. If it must be used, it must be under the supervision of a specially assigned person, or under the guidance of a responsible person (such as a dermatologist).
5. If you suffer from any photosensitivity disorder (photosensitivity), do not use this product. If you suffer from photosensitivity disorders, you may cause serious skin reactions.
6. Do not use on the same day as skincare products containing photosensitive ingredients.
7. Pay attention to protecting your eyes during use, wear eye protectors correctly, and it is recommended to close your eyes during use. If you have a genetic disease with your eyes, do not use a biological mask.
8. If there is any discomfort during the use of this product, please stop using it immediately.
9. This product is not waterproof. Do not immerse the equipment in water or rinse it with water. If you think it has come into contact with water, do not use it.
10. Do not use it for infants or children, and store it out of reach of children.
11. If any parts are damaged, please do not use this product.
12. This product is only for daily use. Do not disassemble, modify, or repair it without authorization.
13. It is recommended to use this product 3-4 times a week, with a single-use time of 10-20 minutes.

